


Disposable Collection Chamber Cartridge

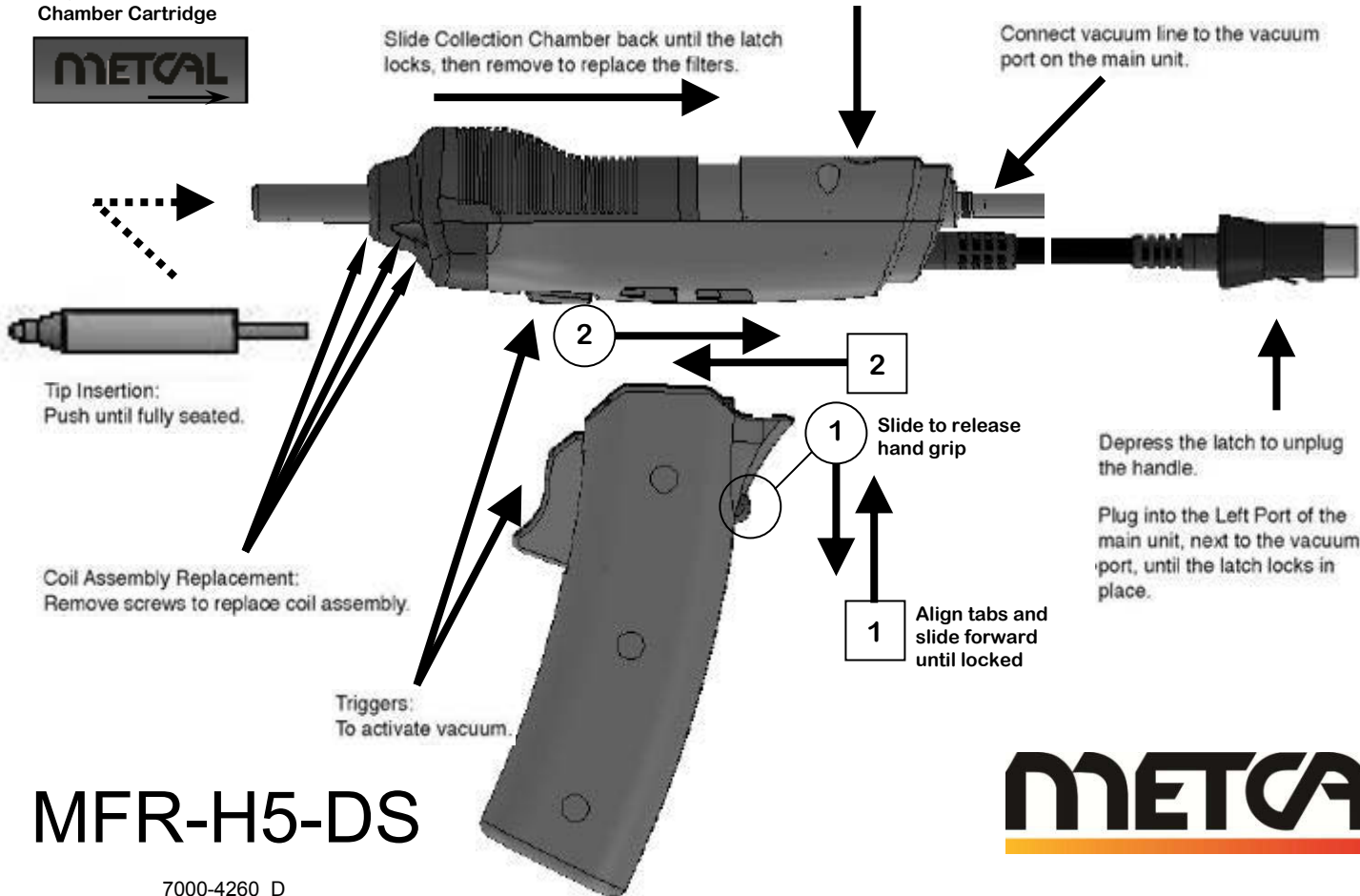

Slide Collection Chamber back until the latch locks, then remove to replace the filters.

Collection Chamber Latch:
Push to lock collection chamber in place.

Connect vacuum line to the vacuum port on the main unit.

Tip Insertion:
Push until fully seated.

Coil Assembly Replacement:
Remove screws to replace coil assembly.

Triggers:
To activate vacuum.

2

1 Slide to release hand grip

1

Align tabs and slide forward until locked

Depress the latch to unplug the handle.

Plug into the Left Port of the main unit, next to the vacuum port, until the latch locks in place.

MFR-H5-DS

7000-4260_D

